

O-Rejuv Facial Device

Lift, Contour, Detox Home Care System

Device and USB Charger

The O-Rejuv is the ultimate homecare system designed to lift, firm, tone and detoxify the skin. This compact, non-invasive, multifunctional device works by combining thermal massage, low level EMS (Electrical Muscle Stimulation) Technology and LED Therapy to stimulate microcirculation to oxygenate the skin for healthier blood flow and nutrient delivery.

Four Treatment Modes:

- **Mode 1: Collagen & Elastin** - Warming Red LED combined with positive microcurrents enhance cellular activity, increase product infusion, support collagen production, reduce fine lines, inflammation and puffiness, restore skin texture and tone, brightens and rejuvenates.
- **Mode 2: Clear & Purify** - Blue LED combined with negative microcurrents provide anti-bacterial benefits, encourage sebum regulation, soothes, detoxifies, clears, calms, improves circulation, brightens, rejuvenates and reduces inflammation and puffiness.
- **Mode 3: Lift & Firm** - Flashing Red LED with low level EMS helps to stimulate collagen and elastin production to firm, lift and tone, increase muscle stimulation, facial contouring, reduce fine lines and improves skin texture and complexion. *Some patients may experience a tingling sensation due to the EMS technology. Adding mode 4 vibration to this setting may decrease this sensation. No thermal energy is emitted during mode 3.*
- **Mode 4: Tone & Plump** - Additional vibration mode can be chosen during any mode, double-click the power button to activate the vibration massage. This will intensify each mode with additional massage to increase microcirculation to oxygenate and rejuvenate the skin.

The O-Rejuv device should not be used by: Children, pregnant women, or people: subject to seizures, with cancers/tumours, with pacemakers, electronic or metal implants (including dental). Seek doctors' clearance for any medical concerns.

A\$199.00 In Stock SKU 00R34

- 1 +

ADD TO CART

FIND A CLINIC



How to use

For optimum results use morning and night for five minutes with your

[ADD TO CART](#)[FIND A CLINIC](#)

How to use

For optimum results use morning and night for five minutes with your favourite O Cosmedics Oil or water-based serum. Perfect to use on the face, neck and décolletage. Care should be taken around the eye area.

Always test this mode on your arm to ensure comfortable heating before applying to the face.

To clean, use a soft, lint-free cloth dampened with water. Ensure the device is off and not plugged into the charger before cleaning. Do not immerse your O-Rejuv device in water at any time.

Best for

For all skins. The O-Rejuv device should not be used by: Children, pregnant women, or people: subject to seizures, with cancers/tumours, with pacemakers, electronic or metal implants (including dental). Seek doctors' clearance for any medical concerns.




How to Build Your Skin Health Tool Kit

Gone are the days where skin care just meant products... As technology has advanced and the world has gotten more and more skin-smart, the rocketing trend of tools and devices has continued to grow. While your homecare tool-kit will never replace the professional strength actives, treatments and results found in an in-clinic treatment, it can work as a brilliant way to enhance and maintain skin health results in between regular treatments and to compliment prescribed and active O COSMEDICS at-home skin routines.


[READ MORE](#)

PAIRS WELL WITH

REAL RESULTS

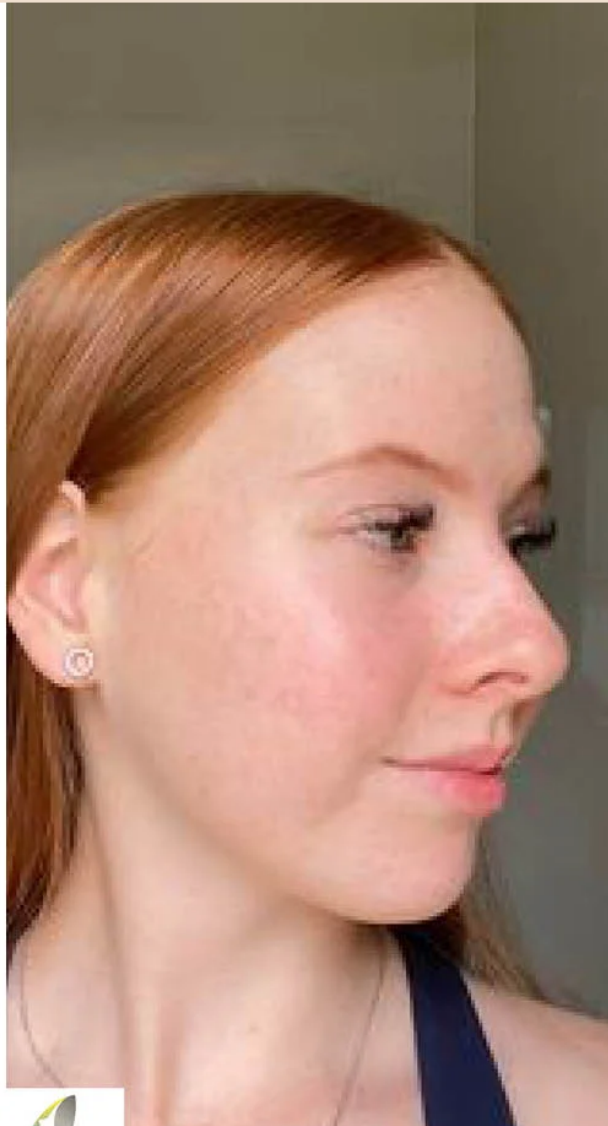

Document title: O-Rejuv - At Home LED Light Therapy Device - O Cosmedics

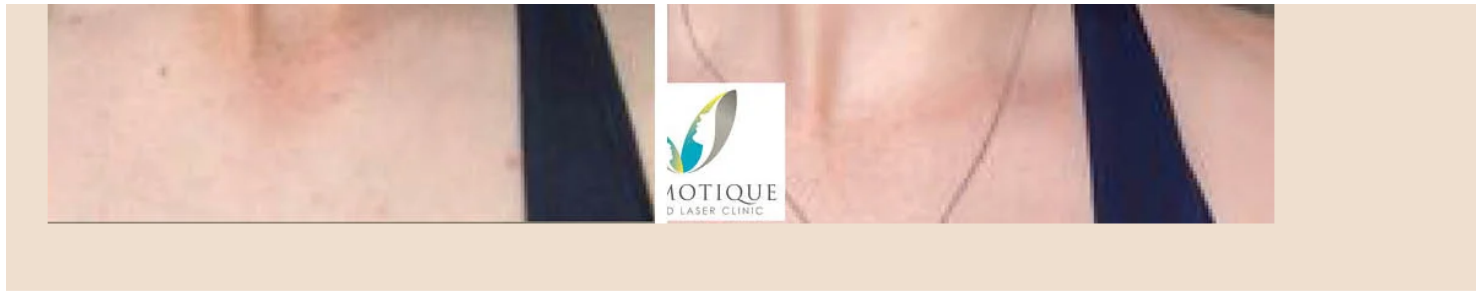

FAQ

What is the O-Rejuv Facial Device?

Is the O-Rejuv Facial Device suitable for all skin types?

How should I incorporate the skin tightening device into my skincare routine?

How does the O-Rejuv Facial Device work?

How often should I use an at-home LED light therapy device?

ABOUT

[ABOUT US](#)

[INGREDIENTS](#)

[O PRO DERMAL PLANNING](#)

[O PRO DERMAL ACTIVE](#)

SHOP

[RETINOL SERUMS](#)

[AHAS AND BHA](#)

[ANTI AGING SKINCARE](#)

[VITAMIN C SKINCARE](#)

[NIACINAMIDE SKINCARE](#)

CUSTOMER CARE

[SHIPPING POLICY](#)

[RETURNS](#)

[TERMS + CONDITIONS](#)

[PRIVACY POLICY](#)

CONTACT US

[CONTACT US](#)

[FIND A CLINIC](#)


[BECOME A PARTNER](#)

Sign up for regular updates & special offers.

NAME

ENTER YOUR EMAIL:

JOIN NOW

 © 2022 O COSMEDICS

JOIN THE O MOVEMENT

